-  EN Molecular Detection Assay 2 - *Listeria*
-  FR Kit de détection moléculaire *Listeria* sp version 2
-  DE Molekulare Detektion 2 – *Listerien* Nachweis
-  IT Analisi molecolare per il rilevamento della *Listeria*
-  ES Ensayo de Detección Molecular 2 para *Listeria*
-  NL Moleculair detectie assay 2 *Listeria*
-  SV Molecular Detection Assay 2 - *Listeria*
-  DA Molekylær Detektions Analyse 2 *Listeria*
-  NO Molekylær deteksjonstest 2 for *Listeria*
-  FI Molekyläärinen testisetti 2 - *Listeria*
-  PT Ensaio para Detecção Molecular de *Listeria* 2
-  EL Δοκιμασία Μοριακής Ανίχνευσης *Listeria* 2
-  PL Molekularny test do wykrywania 2 – *Listeria*
-  RU Молекулярный анализ 2 - *Listeria*
-  TR Moleküler Tayin Testi *Listeri-2*
-  JA 分子検出アッセイ2 - リステリア属菌
-  ZH 李斯特菌分子检测分析 2 - 李斯特菌 (Simplified)
-  ZH 分子檢測套組 2 - 李斯特菌 (Traditional)
-  TH ชุดทดสอบเชื้อลิสทีเรียระดับโนเบกุล 2
-  KO 분자 검출 키트2 - 리스테리아용
-  ID Deteksi Molekuler untuk Pengujian *Listeria* 2

Product Instructions

Molecular Detection Assay 2 - *Listeria*

Product Description and Intended Use

Neogen® Molecular Detection Assay 2 - *Listeria* is used with the Neogen® Molecular Detection System for the rapid and specific detection of *Listeria* species in enriched food and environmental samples.

The Neogen® Molecular Detection Assays use loop-mediated isothermal amplification to rapidly amplify nucleic acid sequences with high specificity and sensitivity, combined with bioluminescence to detect the amplification. Presumptive positive results are reported in real-time while negative results are displayed after the assay is completed. Presumptive positive results should be confirmed using your preferred method or as specified by local regulations^(1, 2, 3).

The Neogen Molecular Detection Assay 2 - *Listeria* is intended for use in a laboratory environment by professionals trained in laboratory techniques. Neogen has not documented the use of this product in industries other than food or beverage. For example, Neogen has not documented this product for testing water, pharmaceutical, cosmetics, clinical or veterinary samples. The Neogen Molecular Detection Assay 2 - *Listeria* has not been evaluated with all possible testing protocols or with all possible strains of bacteria.

As with all test methods, the source, formulation and quality of enrichment medium can influence the results. Factors such as sampling methods, testing protocols, sample preparation, handling, and laboratory technique may also influence results. Neogen recommends evaluation of the method including enrichment medium, in the user's environment using a sufficient number of samples with particular foods and/or environmental samples and microbial challenges to ensure that the method meets the user's criteria. Neogen has evaluated the Neogen Molecular Detection Assay 2 - *Listeria* with Demi-Fraser Broth and Fraser Broth containing Ferric Ammonium Citrate. A typical formulation of these media follows below.

Demi-Fraser Broth Base Typical Formula	(g/L)	Fraser Broth Base Typical Formula	(g/L)
Sodium Chloride	20 g	Sodium Chloride	20 g
Sodium Phosphate, dibasic, anhydrous*	9.6 g	Sodium Phosphate, dibasic, anhydrous*	9.6 g
Beef Extract	5.0 g	Beef Extract	5.0 g
Pancreatic Digest of Casein	5.0 g	Pancreatic Digest of Casein	5.0 g
Peptic Digest of Animal Tissue	5.0 g	Peptic Digest of Animal Tissue	5.0 g
Yeast Extract	5.0 g	Yeast Extract	5.0 g
Lithium Chloride	3.0 g	Lithium Chloride	3.0 g
Potassium Phosphate, monobasic	1.35 g	Potassium Phosphate, monobasic	1.35 g
Esculin	1.0 g	Esculin	1.0 g
Acriflavin HCl	0.0125 g	Acriflavin HCl	0.025 g
Nalidixic Acid	0.01 g	Nalidixic Acid	0.02 g
* Substitute: Sodium Phosphate, dibasic, dihydrate	12.0 g		

Fraser Broth Supplement

(Ingredients per 10 mL vial. One vial is added to one liter of basal medium.)

Ferric Ammonium Citrate 0.5 g/10 mL

Final pH 7.2 ± 0.2 at 25°C

The Neogen® Molecular Detection Instrument is intended for use with samples that have undergone heat treatment during the assay lysis step, which is designed to destroy organisms present in the sample. Samples that have not been properly heat treated during the assay lysis step may be considered a potential biohazard and should NOT be inserted into the Neogen Molecular Detection Instrument.

Neogen Food Safety is certified to ISO (International Organization for Standardization) 9001 for design and manufacturing.

The Neogen Molecular Detection Assay 2 - *Listeria* test kit contains 96 tests, described in Table 1.

**Table 1.** Neogen Molecular Detection Assay Kit Components

Item	Identification	Quantity	Contents	Comments
Neogen® Lysis Solution (LS)	Pink solution in clear tubes	96 (12 strips of 8 tubes)	580 µL of LS per tube	Racked and ready to use
Neogen® Molecular Detection Assay 2 <i>Listeria</i> Reagent Tubes	Blue tubes	96 (12 strips of 8 tubes)	Lyophilized specific amplification and detection mix	Ready to use
Extra caps	Blue caps	96 (12 strips of 8 caps)		Ready to use
Neogen® Reagent Control (RC)	Clear flip-top tubes	16 (2 pouches of 8 individual tubes)	Lyophilized control DNA, amplification and detection mix	Ready to use

The Negative Control, not provided in the kit, is sterile enrichment medium, e.g., Demi-Fraser Broth. Do not use water as a Negative Control.

Safety

The user should read, understand and follow all safety information in the instructions for the Neogen Molecular Detection System and the Neogen Molecular Detection Assay 2 - *Listeria*. Retain the safety instructions for future reference.

⚠ WARNING: Indicates a hazardous situation, which, if not avoided, could result in death or serious injury and/or property damage.

NOTICE: Indicates a potentially hazardous situation which, if not avoided, could result in property damage.

⚠ WARNING

Do not use the Neogen Molecular Detection Assay 2 - *Listeria* in the diagnosis of conditions in humans or animals.

The Neogen Molecular Detection Assay 2 - *Listeria* method may generate *Listeria monocytogenes* to levels sufficient to cause stillbirths and fatalities in pregnant women and the immunocompromised, if exposed.

The user must train its personnel in current proper testing techniques: for example, Good Laboratory Practices, ISO/IEC 17025⁽⁴⁾, or ISO 7218⁽⁵⁾.

To reduce the risks associated with a false-negative result leading to the release of contaminated product:

- Follow the protocol and perform the tests exactly as stated in the product instructions.
- Store the Neogen Molecular Detection Assay 2 - *Listeria* as indicated on the package and in the product instructions.
- Always use the Neogen Molecular Detection Assay 2 - *Listeria* by the expiration date.
- Use the Neogen Molecular Detection Assay 2 - *Listeria* for food and environmental samples that have been validated internally or by a third party.
- Use the Neogen Molecular Detection Assay 2 - *Listeria* only for surfaces, sanitizers, protocols and bacterial strains that have been validated internally or by a third party.
- For an environmental sample containing Neutralizing Buffer with aryl sulfonate complex, perform a 1:2 dilution before testing (1 part sample into 1 part sterile enrichment broth). Another option is to transfer 10 µL of the NB enrichment into the Neogen Lysis Solution tubes. Neogen® Sample Collection products which include Neutralizing Buffer with aryl sulfonate complex: RS96010NB, RS9604NB, SSL10NB, SSL10NB2G, HS10NB, HS10NB2G and HS2410NB2G.

To reduce the risks associated with exposure to chemicals and biohazards:

- It is strongly recommended that female laboratory staff be informed of the risk to a developing fetus resulting from infection of the mother through exposure to *Listeria monocytogenes*.
- Perform pathogen testing in a properly equipped laboratory under the control of trained personnel. Incubated enrichment media and equipment or surfaces that have come into contact with incubated enrichment media may contain pathogens at levels sufficient to cause risk to human health.
- Always follow standard laboratory safety practices, including wearing appropriate protective apparel and eye protection while handling reagents and contaminated samples.
- Avoid contact with the contents of the enrichment media and reagent tubes after amplification.

- Dispose of enriched samples according to current industry standards.
- Samples that have not been properly heat treated during the assay lysis step may be considered a potential biohazard and should NOT be inserted into the Neogen Molecular Detection Instrument.

To reduce the risks associated with cross-contamination while preparing the assay:

- Always wear gloves (to protect the user and prevent introduction of nucleases).

To reduce the risks associated with environmental contamination:

- Follow current industry standards for disposal of contaminated waste.

To reduce the risks associated with exposure to hot liquids:

- Do not exceed the recommended temperature setting on heater.
- Do not exceed the recommended heating time.
- Use an appropriate, calibrated thermometer to verify the Neogen® Molecular Detection Heat Block Insert temperature (e.g., a partial immersion thermometer or digital thermocouple thermometer, not a total immersion thermometer.) The thermometer must be placed in the designated location in the Neogen Molecular Detection Heat Block Insert.

NOTICE**To reduce the risks associated with cross-contamination while preparing the assay:**

- Change gloves prior to reagent pellet hydration.
- Use of sterile, aerosol barrier (filtered), molecular biology grade pipette tips is recommended.
- Use a new pipette tip for each sample transfer.
- Use Good Laboratory Practices to transfer the sample from the enrichment to the lysis tube. To avoid pipettor contamination, the user may choose to add an intermediate transfer step. For example, the user can transfer each enriched sample into a sterile tube.
- Use a molecular biology workstation containing germicidal lamp where available. Periodically decontaminate laboratory benches and equipment (pipettes, cap/decap tools, etc.) with a 1- 5% (v:v in water) household bleach solution or DNA removal solution.

To reduce the risks associated with a false-positive result:

- Never open reagent tubes post amplification.
- Always dispose of the contaminated tubes by soaking in a 1-5% (v:v in water) household bleach solution or equivalent solution for 1 hour and away from the assay preparation area.
- Never autoclave reagent tubes post amplification.

Consult the Safety Data Sheet for additional information and local regulations for disposal.

If you have questions about specific applications or procedures, please visit our website at www.neogen.com or contact your local Neogen representative or distributor.

User Responsibility

Users are responsible for familiarizing themselves with product instructions and information. Visit our website at www.neogen.com, or contact your local Neogen representative or distributor for more information.

When selecting a test method, it is important to recognize that external factors such as sampling methods, testing protocols, sample preparation, handling, and laboratory technique may influence results.

It is the user's responsibility in selecting any test method or product to evaluate a sufficient number of samples with the appropriate matrices and microbial challenges to satisfy the user that the chosen test method meets the user's criteria.

It is also the user's responsibility to determine that any test methods and results meet its customers' and suppliers' requirements.

As with any test method, results obtained from use of any Neogen Food Safety product do not constitute a guarantee of the quality of the matrices or processes tested.

To help customers evaluate the method for various food matrices, Neogen has developed the Neogen® Molecular Detection Matrix Control kit. When needed, use the Matrix Control (MC) to determine if the matrix has the ability to impact the Neogen Molecular Detection Assay 2 - *Listeria* results. Test several samples representative of the matrix, i.e., samples obtained from different origin, during any validation period when adopting the Neogen method or when testing new or unknown matrices or matrices that have undergone raw material or process changes.

A matrix can be defined as a type of product with intrinsic properties such as composition and process. Differences between matrices may be as simple as the effects caused by differences in their processing or presentation for example, raw vs. pasteurized; fresh vs. dried, etc.

Limitation of Warranties / Limited Remedy

EXCEPT AS EXPRESSLY STATED IN A LIMITED WARRANTY SECTION OF INDIVIDUAL PRODUCT PACKAGING, NEOGEN DISCLAIMS ALL EXPRESS AND IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO, ANY WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR USE. If any Neogen Food Safety Product is defective, Neogen or its authorized distributor will, at its option, replace or refund the purchase price of the product. These are your exclusive remedies. Please contact your Neogen representative or authorized Neogen distributor for any further questions.

Limitation of Neogen Liability

NEOGEN WILL NOT BE LIABLE FOR ANY LOSS OR DAMAGES, WHETHER DIRECT, INDIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES, INCLUDING BUT NOT LIMITED TO LOST PROFITS. In no event shall Neogen's liability under any legal theory exceed the purchase price of the product alleged to be defective.

Storage and Disposal

Store the Neogen Molecular Detection Assay 2 - *Listeria* at 2-8°C. Do not freeze. Keep kit away from light during storage. After opening the kit, check that the foil pouch is undamaged. If the pouch is damaged, do not use. After opening, unused reagent tubes should always be stored in the re-sealable pouch with the desiccant inside to maintain stability of the lyophilized reagents. Store resealed pouches at 2-8°C for no longer than 60 days.

Do not use Neogen Molecular Detection Assay 2 - *Listeria* past the expiration date. Expiration date and lot number are noted on the outside label of the box. After use, the enrichment medium and the Neogen Molecular Detection Assay 2 - *Listeria* tubes can potentially contain pathogenic materials. When testing is complete, follow current industry standards for the disposal of contaminated waste. Consult the Safety Data Sheet for additional information and local regulations for disposal.

Instructions for Use

Follow all instructions carefully. Failure to do so may lead to inaccurate results.

Periodically decontaminate laboratory benches and equipment (pipettes, cap/decap tools, etc.) with a 1- 5% (v:v in water) household bleach solution or DNA removal solution.

The user should complete the Neogen Molecular Detection System operator qualification training, as described in the "Installation Qualification (IQ) / Operational Qualification (OQ) Protocols and Instructions for Neogen Molecular Detection System" document⁽⁶⁾.

See Section "Specific Instructions for validated methods" for specific requirements:

Table 3 for enrichment protocols according to AOAC® Official Method of AnalysisSM 2016.07 and AOAC® Performance Tested MethodSM Certificate #111501.

Table 4 for enrichment protocols according to NF Validation certificate 3M 01/14-05/16.

Sample Enrichment

Tables 2, 3 or 4 present guidance for the enrichment of food and environmental samples. It is the user's responsibility to validate alternative sampling protocols or dilution ratios to ensure this test method meets the user's criteria.

Foods

1. Allow the Demi-Fraser Broth enrichment medium (includes ferric ammonium citrate) to equilibrate to ambient laboratory temperature.
2. Aseptically combine the enrichment medium and sample according to Tables 2, 3 or 4. For all meat and highly particulate samples, the use of filter bags is recommended.
3. Homogenize thoroughly by blending, stomaching, hand mixing for 2 ± 0.2 minutes. Incubate at $37 \pm 1^\circ\text{C}$ according to Tables 2, 3 or 4.
4. For raw dairy products (see Tables 2 or 4), transfer 0.1 mL of the primary enrichment into 10 mL of Fraser Broth. Incubate at $37 \pm 1^\circ\text{C}$ for 20-24 hours.

Environmental Samples

Sample collection devices can be a sponge hydrated with a neutralizing solution to inactivate the effects of the sanitizers.

Neogen recommends the use of a biocide-free cellulose sponge. Neutralizing solution can be Dey-Engley (D/E) Neutralizing Broth or Letheen Broth. It is recommended to sanitize the area after sampling.

WARNING: Should you select to use Neutralizing Buffer (NB) that contains aryl sulfonate complex as the hydrating solution for the sponge, it is required to perform a 1:2 dilution (1 part sample into 1 part sterile enrichment broth) of the enriched environmental sample before testing in order to reduce the risks associated with a false-negative result leading to the release of contaminated product. Another option is to transfer 10 µL of the neutralizing buffer enrichment into the Neogen Lysis Solution tubes.

The recommended size of the sampling area to verify the presence or absence of the pathogen on the surface is at least 100 cm² (10 cm x 10 cm or 4"x4"). When sampling with a sponge, cover the entire area going in two directions (left to right then up and down) or collect environmental samples following your current sampling protocol or according to the FDA BAM⁽¹⁾, USDA FSIS MLG⁽²⁾ or ISO 18593⁽⁷⁾ guidelines.

1. Allow the Demi-Fraser Broth enrichment medium (includes ferric ammonium citrate) to equilibrate to ambient laboratory temperature.
2. Aseptically combine the enrichment medium and sample according to Tables 2, 3 or 4.
3. Homogenize thoroughly by blending, stomaching, hand mixing or vortexing for 2 ± 0.2 minutes. Incubate at 37 ± 1°C for 24-30 hours according to Tables 2, 3 or 4.

Table 2: General enrichment protocols at 37 ± 1°C using Demi-Fraser Broth^(a) and Fraser Broth^(a) as needed.

Sample Matrix	Sample Size	Enrichment Broth Volume (mL)	Enrichment Temperature (± 1°C)	Enrichment Time (hours)
Heat-processed, cooked, cured meats, poultry, seafood and fish				
Heat-processed / pasteurized dairy products	25 g	225	37	24-30
Produce and vegetables				
Multi-component foods				
Environmental samples	1 sponge	100 or 225	37	24-30
	1 swab	10	37	24-30
Raw meat, poultry, seafood, fish	25 g	475	37	28-32

Sample Matrix	Primary Enrichment (Demi-Fraser Broth) ^(a)				Secondary Enrichment (Fraser Broth) ^(a)			Sample Analysis Volume ^(b)
	Sample Size	Enrichment Broth Volume (mL)	Enrichment Temperature ($\pm 1^{\circ}\text{C}$)	Enrichment Time (hours)	Sample Size	Enrichment Temperature ($\pm 1^{\circ}\text{C}$)	Enrichment Time (hours)	
Raw dairy products	25 g	225	37	20-24	Transfer 0.1 mL into 10 mL Fraser Broth	37	20-24	10 μL

(a) Demi-Fraser and Fraser Broth should always be supplemented with Fraser Broth Supplement (ferric ammonium citrate) during primary or secondary enrichment.

(b) Volume of enriched sample transferred to the Neogen Lysis Solution tubes. Refer to Step 4.6 of Lysis section.

Specific Instructions for Validated Methods

AOAC® Official MethodSM 2016.07

AOAC® Performance Tested MethodSM #111501

In AOAC Research Institute OMASM and PTMSM programs, the Neogen Molecular Detection Assay 2 – *Listeria* was found to be an effective method for the detection of *Listeria* species. The matrices tested in the study are shown in Table 3.

Table 3. Enrichment protocols using Demi-Fraser Broth^(a) at $37 \pm 1^{\circ}\text{C}$ according to AOAC OMASM 2016.07 and AOAC PTMSM Certificate #111501.

Sample Matrix		Sample Size	Enrichment Broth Volume (mL)	Enrichment Time (hours)
Beef hot dogs, queso fresco, vanilla ice cream, 4% milk fat cottage cheese, 3% chocolate whole milk, romaine lettuce, bagged raw spinach, cold smoked salmon	25 g		225	24-30
Raw chicken	25 g		475	28-32
Deli turkey	125 g		1125	24-30
Cantaloupe ^(b)	Whole melon	Enough volume to allow melon to float		26-30
Environmental samples:	Stainless steel	1 sponge	225	24-30
	Sealed concrete	1 sponge	100	24-30
	Plastic ^(c)	1 swab	10	24-30

All samples for the AOAC validation were homogenized by stomaching unless otherwise noted.

(a) Demi-Fraser and Fraser Broth should always be supplemented with Fraser Broth Supplement (ferric ammonium citrate) during primary or secondary enrichment.

(b) Homogenize sample by hand mixing.

(c) Homogenize sample by vortexing.

NF Validation by AFNOR Certification

3M 01/14-05/16

ALTERNATIVE ANALYTICAL METHODS FOR AGRIBUSINESS

<http://nf-validation.afnor.org/en>

For more information about end of validity, please refer to NF VALIDATION certificate available on the website mentioned above.

NF VALIDATION certified method in compliance with ISO 16140-2⁽⁸⁾ in comparison to ISO 11290-1⁽³⁾

Scope of the validation: All human food and environmental samples (excluding primary production samples)

Sample preparation: Samples should be prepared according to ISO 11290-1⁽³⁾ and ISO 6887⁽⁹⁾

Software version: See certificate

Table 4. Enrichment protocols according to NF VALIDATION certified method 3M 01/14-05/16 at 37 ± 1 °C using Demi-Fraser Broth^(a) and Fraser Broth^(a) as needed.

General Protocol	Sample Size	Enrichment Broth Volume (mL)	Enrichment Temperature (± 1°C)	Enrichment Time (hours)	Sample Analysis Volume ^(b)	Accepted Interruption Point
All food samples (excepted raw meats, raw seafood and raw dairy products)	25 g	225	37	24-30	20 µL	<ul style="list-style-type: none"> • Demi-Fraser Broth up to 72 hours • Lysate at -20°C • Lysate at 4°C up to 72 hours
Environmental samples	25 g, 1 swab, or 1 wipe	225	37	24-30	20 µL	<ul style="list-style-type: none"> • Lysate at -20°C
Specific Protocol 1	Sample Size	Enrichment Broth Volume (mL)	Enrichment Temperature (± 1°C)	Enrichment Time (hours)	Sample Analysis Volume ^(b)	Accepted Interruption Point
Raw meats and raw seafood	25 g	475	37	28-32	20 µL	<ul style="list-style-type: none"> • Demi-Fraser Broth up to 72 hours • Lysate at -20°C • Lysate at 4°C up to 72 hours

Specific Protocol 2	Primary Enrichment (Demi-Fraser Broth) ^(a)				Secondary Enrichment (Fraser Broth) ^(a)				Accepted Interruption Point
	Sample Size	Enrichment Broth Volume (mL)	Enrichment Temperature ($\pm 1^{\circ}\text{C}$)	Enrichment Time (hours)	Sample Size	Enrichment Temperature ($\pm 1^{\circ}\text{C}$)	Enrichment Time (hours)	Sample Analysis Volume ^(b)	
Raw dairy products	25 g	225	37	20-24	Transfer 0.1 mL into 10 mL Fraser Broth	37	20-24	10 μL	<ul style="list-style-type: none"> Fraser Broth up to 72 hours Lysate at -20°C

(a) Demi-Fraser and Fraser Broth should always be supplemented with Fraser Broth Supplement (ferric ammonium citrate) during primary or secondary enrichment.

(b) Volume of enriched sample transferred to the Neogen Lysis Solution tubes. Refer to Step 4.6 of Lysis section.

Note: Samples larger than 25 g have not been tested in the NF validation study.

Preparation of the Neogen® Molecular Detection Speed Loader Tray

1. Wet a cloth with a 1-5% (v:v in water) household bleach solution or equivalent solution and wipe the Neogen® Molecular Detection Speed Loader Tray.
2. Rinse the Neogen Molecular Detection Speed Loader Tray with water.
3. Use a disposable towel to wipe the Neogen Molecular Detection Speed Loader Tray dry.
4. Ensure the Neogen Molecular Detection Speed Loader Tray is dry before use.

Preparation of the Neogen® Molecular Detection Chill Block Insert

Place the Neogen® Molecular Detection Chill Block directly on the laboratory bench. Use the chill block at ambient laboratory temperature (20-25°C).

Preparation of the Neogen® Molecular Detection Heat Block Insert

Place the Neogen® Molecular Detection Heat Block Insert in a dry double block heater unit. Turn on the dry block heater unit and set the temperature to allow the Neogen Molecular Detection Heat Block Insert to reach and maintain a temperature of 100 $\pm 1^{\circ}\text{C}$.

Note: Depending on the heater unit, allow approximately 30 minutes for the Neogen Molecular Detection Heat Block Insert to reach temperature. Using an appropriate, calibrated thermometer (e.g., a partial immersion thermometer or digital thermocouple thermometer, not a total immersion thermometer) placed in the designated location, verify that the Neogen Molecular Detection Heat Block Insert is at 100 $\pm 1^{\circ}\text{C}$.


Preparation of the Neogen® Molecular Detection Instrument

1. Launch the Neogen® Molecular Detection System Software and log in. Contact your Neogen Food Safety representative to ensure you have the most updated version of the software.
2. Turn on the Neogen Molecular Detection Instrument.
3. Create or edit a run with data for each sample. Refer to the Neogen Molecular Detection System User Manual for details.

Note: The Neogen Molecular Detection Instrument must reach Ready State before inserting the Neogen Molecular Detection Speed Loader Tray with reaction tubes. This heating step takes approximately 20 minutes and is indicated by an ORANGE light on the instrument's status bar. When the instrument is ready to start a run, the status bar will turn GREEN.


Lysis

Remove the bottom of Neogen Lysis Solution Rack with a screwdriver or spatula before placing in the Neogen Molecular Detection Heat Block Insert.

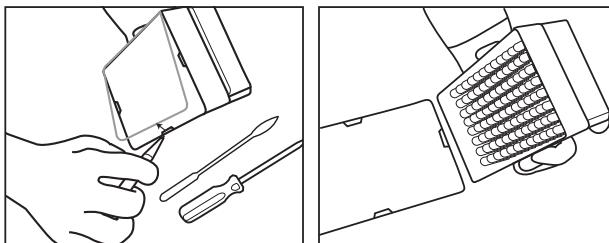1. Allow the Neogen Lysis Solution tubes to warm up by setting the rack at room temperature (20-25°C) overnight (16-18 hours). Alternatives to equilibrate the Neogen Lysis Solution tubes to room temperature are to set the Neogen Lysis Solution tubes on the laboratory bench for at least 2 hours, incubate the Neogen Lysis Solution tubes in a 37 ± 1°C incubator for 1 hour or place them in a dry double block heater for 30 seconds at 100 ± 1°C.
2. Invert the capped tubes to mix. Proceed to next step within 4 hours.
3. Remove the enrichment broth from the incubator.
4. One Neogen Lysis Solution tube is required for each sample and the Negative Control (NC) sample (sterile enrichment medium).
 - 4.1 Neogen Lysis Solution tube strips can be cut to desired Neogen Lysis Solution tube number. Select the number of individual Neogen Lysis Solution tubes or 8-tube strips needed. Place the Neogen Lysis Solution tubes in an empty rack.
 - 4.2 To avoid cross-contamination, decap one Neogen Lysis Solution tube strip at a time and use a new pipette tip for each transfer step.
 - 4.3 Transfer enriched sample to Neogen Lysis Solution tubes as described below:

Transfer each enriched sample into an individual Neogen Lysis Solution tube **first**. Transfer the NC **last**.
- 4.4 Use the Neogen® Molecular Detection Cap/Decap Tool-Lysis to decap one Neogen Lysis Solution tube strip - one strip at a time.
- 4.5 Discard the Neogen Lysis Solution tube cap – if lysate will be retained for retest, place the caps into a clean container for re-application after lysis.
 - 4.5.1 For processing of retained lysate, see Appendix A.
- 4.6 Transfer 20 µL of sample into a Neogen Lysis Solution tube unless otherwise indicated in protocols from Tables 2, 3 or 4 (e.g., raw dairy products or when using environmental samples with neutralizing buffer use 10 µL).

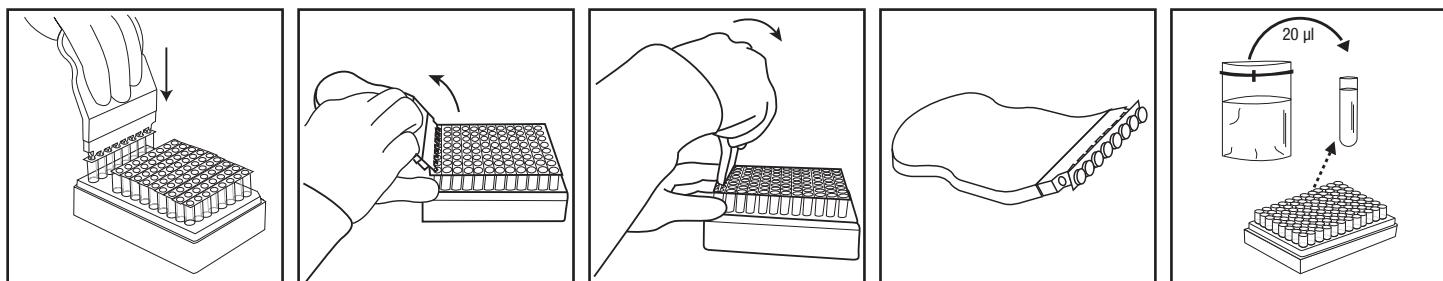
5. Repeat step 4.4 to 4.6 until each individual sample has been added to a corresponding Neogen Lysis Solution tube in the strip.
6. When all samples have been transferred, transfer 20 µL of NC (sterile enrichment medium, e.g., Demi-Fraser Broth) into a Neogen Lysis Solution tube. Do not use water as a NC.
7. Verify that the temperature of the Neogen Molecular Detection Heat Block Insert is at 100 ± 1°C.
8. Place the uncovered rack of Neogen Lysis Solution tubes in the Neogen Molecular Detection Heat Block Insert and heat for 15 ± 1 minutes. During heating, the Neogen Lysis Solution will change from pink (cool) to yellow (hot).
 - 8.1 Samples that have not been properly heat treated during the assay lysis step may be considered a potential biohazard and should NOT be inserted into the Neogen Molecular Detection Instrument.

9. Remove the uncovered rack of Neogen Lysis Solution tubes from the heating block and allow to cool in the Neogen Molecular Detection Chill Block Insert at least 5 minutes and a maximum of 10 minutes. The Neogen Molecular Detection Chill Block Insert used at ambient temperature should sit directly on the laboratory bench. When cool, the lysis solution will revert to a pink color.

10. Remove the rack of Neogen Lysis Solution tubes from the Neogen Molecular Detection Chill Block Insert.

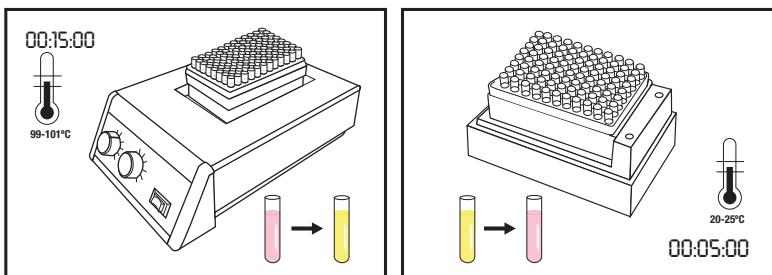
Amplification

1. One Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tube is required for each sample and the NC.
 - 1.1 Reagent tube strips can be cut to desired tube number. Select the number of individual Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tubes or 8-tube strips needed.
 - 1.2 Place Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tubes in an empty rack.
 - 1.3 Avoid disturbing the reagent pellets from the bottom of the tubes.
2. Select one Neogen Reagent Control (RC) tube and place in rack.
3. To avoid cross-contamination, decap one Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tube strip at a time and use a new pipette tip for each transfer step.
4. Transfer lysate to Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tubes and Neogen Reagent Control tube as described below:

Transfer each sample lysate into individual Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tubes **first** followed by the NC. Hydrate the RC tube **last**.

5. Use the Neogen® Molecular Detection Cap/Decap Tool-Reagent to decap the Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tubes –one Reagent tube strip at a time. Discard cap.
 - 5.1 Transfer 20 µL of sample lysate from the upper 1/2 of the liquid (avoid precipitate) in the Neogen Lysis Solution tube into corresponding Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tube. Dispense at an angle to avoid disturbing the pellets. Mix by gently pipetting up and down 5 times.
 - 5.2 Repeat step 5.1 until individual sample lysate has been added to a corresponding Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tube in the strip.
 - 5.3 Cover the Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tubes with the provided extra caps and use the rounded side of the Neogen Molecular Detection Cap/Decap Tool-Reagent to apply pressure in a back and forth motion ensuring that the cap is tightly applied.
 - 5.4 Repeat steps 5.1 and 5.3 as needed, for the number of samples to be tested.
 - 5.5 When all sample lysates have been transferred, repeat steps 5.1 and 5.3 to transfer 20 µL of NC lysate into a Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tube.
 - 5.6 Transfer **20 µL of NC lysate** into a Neogen Reagent Control tube. Dispense at an angle to avoid disturbing the pellets. Mix by gently pipetting up and down 5 times.
6. Load capped tubes into a clean and decontaminated Neogen Molecular Detection Speed Loader Tray. Close and latch the Neogen Molecular Detection Speed Loader Tray lid.

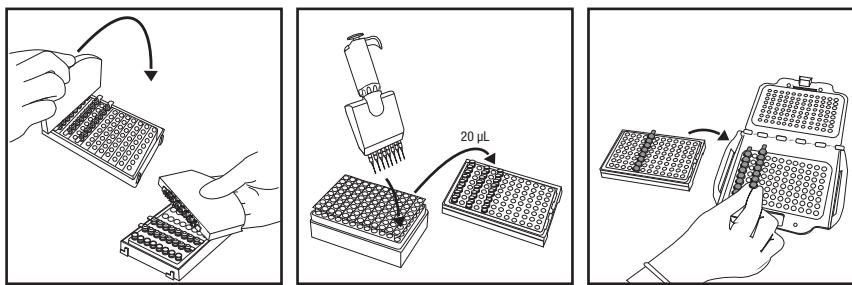7. Review and confirm the configured run in the Neogen Molecular Detection System Software.
8. Click the Start button in the software and select instrument for use. The selected instrument's lid automatically opens.
9. Place the Neogen Molecular Detection Speed Loader Tray into the Neogen Molecular Detection Instrument and close the lid to start the assay. Results are provided within 75 minutes, although positives may be detected sooner.
10. After the assay is complete, remove the Neogen Molecular Detection Speed Loader Tray from the Neogen Molecular Detection Instrument and dispose of the tubes by soaking in a 1-5% (v:v in water) household bleach solution or equivalent solution for 1 hour and away from the assay preparation area.

NOTICE: To minimize the risk of false positives due to cross-contamination, never open Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tubes containing amplified DNA. This includes Neogen Molecular Detection Assay 2 – *Listeria* Reagent Tubes, Neogen Reagent Control and Neogen Molecular Detection Matrix Control tubes. Always dispose of sealed reagent tubes by soaking in a 1-5% (v:v in water) household bleach solution or equivalent solution for 1 hour and away from the assay preparation area.

Results and Interpretation

An algorithm interprets the light output curve resulting from the detection of the nucleic acid amplification. Results are analyzed automatically by the software and are color-coded based on the result. A Positive or Negative result is determined by analysis of a number of unique curve parameters. Presumptive positive results are reported in real-time while Negative and Inspect results will be displayed after the run is completed.

Presumptive positive samples should be confirmed as per the laboratory standard operating procedures or by following the appropriate reference method confirmation^(1, 2, 3), beginning with transfer from the primary enrichment to secondary enrichment broth (if applicable), followed by subsequent plating and confirmation of isolates using appropriate biochemical and serological methods.

NOTE: Even a negative sample will not give a zero reading as the system and Neogen Molecular Detection Assay 2 - *Listeria* amplification reagents have a “background” relative light unit (RLU) reading.

In the rare event of any unusual light output, the algorithm labels this as “Inspect.” Neogen recommends the user to repeat the assay for any Inspect samples. If the result continues to be Inspect, proceed to confirmation test using your preferred method or as specified by local regulations.

Table 5: Result symbols for samples and interpretation as provided by the Neogen Molecular Detection System Software.

Well Type	Well Result Symbol	Result	Interpretation
Sample		Positive	The sample is presumptive positive for the target pathogen.
Sample		Negative	The sample is negative for the target pathogen.
Sample		Inhibited	The sample matrix was inhibitory to the assay. A re-test may be required. Refer to the troubleshooting section and the assay kit Product Instructions for more information.
Sample		Inspect	The presence or absence of the target pathogen was indeterminate. A re-test may be required. Refer to the troubleshooting section and the assay kit Product Instructions for more information.
Sample		Error	No bioluminescence was detected. A re-test may be required. Refer to the troubleshooting section and the assay kit Product Instructions for more information.

Confirmation of Results According to the NF VALIDATION Certified Method

In the context of the NF VALIDATION, all samples identified as positive by the Neogen Molecular Detection Assay 2 - *Listeria* must be confirmed by one of the following tests:

Option 1: Using the ISO 11290-1⁽³⁾ standard starting from Demi-Fraser enrichment.

Option 2: By isolation on PALCAM agar or on chromogenic agar forming part of a NF VALIDATION certified method for detection of the *Listeria* genus. The presence of characteristic colonies is sufficient to confirm the presence of *Listeria* spp.

Option 3: Using nucleic acid probes as described in ISO 7218⁽⁵⁾ standard, performed on isolated colonies, from selective agar (see Options 1 or 2).

Option 4: Using any other method certified NF VALIDATION, the principle of which must be different from Neogen Molecular Detection Assay 2 – *Listeria*. The complete protocol described for this second validated method must be used. All steps prior to the start of confirmation must be common to both methods.

In the event of discordant results (presumptive positive with the Neogen Molecular Detection Assay 2 – *Listeria*, non-confirmed by one of the means described above) the laboratory must follow the necessary steps to ensure the validity of the results obtained.

If you have questions about specific applications or procedures, please visit our website at www.neogen.com or contact your local Neogen representative or distributor.

Appendix A. Protocol Interruption: Storage and re-testing of heat-treated lysates. The storage of lysates at 2 to 8°C for 72 hours does not apply to dairy and environmental samples.

1. To store a heat-treated lysate, re-cap the lysis tube with a clean cap (see “LYSIS”, 4.5)
2. Store at 2 to 8°C for up to 72 hours.
3. Prepare a stored sample for amplification by inverting 2-3 times to mix.
4. Decap the tubes.
5. Place the mixed lysate tubes on Neogen Molecular Detection Heat Block Insert and heat at 100 ± 1°C for 5 ± 1 minutes.
6. Remove the rack of Neogen Lysis Solution tubes from the heating block and allow to cool in the Neogen Molecular Detection Chill Block Insert at least 5 minutes and a maximum of 10 minutes.
7. Continue the protocol at the ‘Amplification’ section detailed above.

References:

1. US Food and Drug Administration Bacteriological Analytical Manual. Chapter 10: Detection of *Listeria monocytogenes* in foods and environmental samples, and enumeration of *Listeria monocytogenes* in foods. October 31, 2017.
2. US Department of Agriculture (USDA) FSIS Microbiology Laboratory Guidebook 8.11. Isolation and identification of *Listeria monocytogenes* from red meat, poultry, ready-to-eat siluriformes (fish) and egg products, and environmental samples. January 2, 2019.
3. NF EN ISO 11290-1. Microbiology of the food chain – Horizontal method for the detection and enumeration of *Listeria monocytogenes* and of *Listeria* spp. – Part 1: Detection method.
4. NF EN ISO/IEC 17025. General requirements for the competence of testing and calibration laboratories.
5. NF EN ISO 7218. Microbiology of food and animal feeding stuffs – General requirements and guidance for microbiological examinations.
6. Neogen Installation Qualification (IQ) / Operational Qualification (OQ) Protocols and Instructions for Neogen Molecular Detection System.
7. NF EN ISO 18593. Microbiology of the food chain – Horizontal methods for surface sampling.
8. NF EN ISO 16140-2 Microbiology of the food chain – Method validation – Part 2: Protocol for the validation of alternative (proprietary) methods against a reference method.
9. NF EN ISO 6887. Microbiology of the food chain – Preparation of test samples, initial suspension and decimal dilutions for microbiological examination.

Explanation of Product Label Symbols

info.neogen.com/symbols

Neogen Food Safety

Neogen Corporation
620 Lesher Place
Lansing, MI 48912 USA
Neogen.com

Neogen Europe Ltd.
The Dairy School
Auchincruive
Ayr, KA6 5HU
Scotland, UK

Neogen Ireland, Ltd.
Bray Business Park, Bray
Co. Wicklow
A98YV29, Ireland

Neogen Corporation
620 Lesher Place Lansing, MI 48912 USA
www.neogen.com

© Neogen Corporation 2024. All rights reserved.
Neogen is a trademark of Neogen Corporation.
All other trademarks are the property of their respective companies.

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