700007660 | MDALRN96

Product Instructions

-  **(EN)** Neogen® Molecular Detection Assay – *Listeria* Right Now™
-  **(FR)** Kit de détection moléculaire Neogen® – *Listeria* Right Now™
-  **(DE)** Neogen® Molekularer Detektions-Assay – *Listeria* Right Now™
-  **(ES)** Neogen®, ensayo de detección molecular: *Listeria* Right Now™
-  **(PT)** Ensaio para Detecção Molecular Neogen® – *Listeria* Right Now™
-  **(JP)** Neogen® 病原菌検出アッセイ Right Now (TM) リステリア属菌 環境モニタリング用
-  **(ZH)** Neogen® 分子检测分析 – *Listeria* Right Now™
-  **(TW)** Neogen® 分子檢測系統 – *Listeria* Right Now™
-  **(TH)** ชุดทดสอบเชื้อจุลินทรีย์ก่อโรคระดับโมเลกุลของ Neogen® – *Listeria* Right Now™
-  **(KO)** Neogen® 분자 검출 분석 – *Listeria* Right Now™

LRN

Listeria Right Now™

Product Instructions

Neogen® Molecular Detection Assay – *Listeria* Right Now™

Product Description and Intended Use

The Neogen Molecular Detection Assay – *Listeria* Right Now (MDALRN96) kit is used with the Neogen® Molecular Detection System for rapid and specific detection of *Listeria* spp. in samples from visibly clean post-cleaning/sanitation and pre-production environmental surfaces without the need for enrichment. It uses the included Neogen Right Now™ Sampler, wet with a sampling solution, to monitor *Listeria* in the environment from food and beverage manufacturing facilities.

The screening kit contains reagent tubes with blue, color-coded lyophilized pellets to detect *Listeria* species strains.

After sample collection using the Right Now Sampler, the contents are lysed to release nucleic acids. Ribosomal RNA is then reverse transcribed into complementary DNA (cDNA). The assay employs loop-mediated isothermal amplification (LAMP) for rapid, highly specific, and sensitive nucleic acid amplification, with bioluminescent detection of the amplified product.

The Neogen Molecular Detection Assay – *Listeria* Right Now is intended for use in a laboratory environment by professionals trained in laboratory techniques. Neogen has not documented the use of this product in industries other than food or beverage. For example, Neogen has not documented this product for testing pharmaceutical, cosmetics, clinical, or veterinary samples. The Neogen Molecular Detection Assay – *Listeria* Right Now has not been evaluated with all possible food process environments or with all possible strains of bacteria. Neogen recommends evaluation of the method in the user's environment using a sufficient number of samples for particular environmental challenges to ensure that the method meets the user's criteria.

Use the paired Neogen® Right Now™ Matrix Control (RNMC96) to demonstrate assay robustness in the user environment when implementing the method or during method verification.

The Neogen® Molecular Detection Instrument is intended for use with samples that have undergone heat treatment during the assay lysis step, which is designed to lyse the cells and release the nucleic acids from the organisms present in the sample. Samples that have not been properly heat treated during the assay lysis step may be considered a potential biohazard and should NOT be inserted into the Neogen Molecular Detection Instrument.

Neogen Food Safety is certified ISO (International Organization for Standardization) 9001 for design and manufacturing.

The Neogen Molecular Detection Assay – *Listeria* Right Now test kit contains enough reagents for 96 tests, as described in Table 1.

Table 1. Kit Components

Item	Identification	Quantity	Contents	Comments
Right Now™ Sampler	Blue capped tubes with polyurethane samplers	96 samplers	Sampler wet with 1.5 mL sampling solution	Ready to use
Right Now™ Lysis Reagent Resuspension Buffer	Clear plastic bottles	2 bottles (Each can reconstitute two lysis buffer reagent bottles)	Resuspension buffer	Ready to use
Right Now™ Lysis Buffer Reagent	Clear glass bottles	4 bottles (Reconstituted, each can fill 28 lysis tubes)	Lysis reagent	Lyophilized
Lysis (LS) tubes	Clear tubes	96 (12 strips of 8 tubes each)	Empty	Racked and ready to use
<i>Listeria</i> Right Now™ Reagent Tubes	Clear tubes with blue pellets	96 (3 pouches containing 4 strips of 8 tubes)	Lyophilized specific amplification and detection mix	Ready to use
Extra caps	Clear caps	96 (12 strips of 8 caps)		Ready to use
Right Now™ Reagent Control (RC)	Clear flip-top tubes	16 (2 pouches of 8 individual tubes)	Lyophilized control nucleic acid, amplification, and detection mix	Ready to use

Safety

The user should read, understand, and follow all safety information in the instructions for the Neogen Molecular Detection System and the Neogen Molecular Detection Assay – *Listeria* Right Now. Retain the safety instructions for future reference.

	WARNING:	Indicates a hazardous situation, which, if not avoided, could result in death or serious injury and/or property damage.
	CAUTION:	Indicates a hazardous situation, which, if not avoided, could result in minor or moderate injury and/or property damage.

WARNING

Do not use the Neogen Molecular Detection Assay – *Listeria* Right Now™ in the diagnosis of conditions in humans or animals.

The user must train personnel in current proper testing techniques, for example: Good Laboratory Practices, ISO/IEC 17025(1), or ISO 7218(2).

The Neogen Molecular Detection Assay – *Listeria* Right Now is not designed for sampling soiled surfaces or for food matrix testing.

To reduce the risks associated with a false-negative result leading to the release of contaminated surfaces:

- Follow the protocol and perform the tests exactly as stated in the product instructions.
- Always use Good Laboratory Practices to transfer samples and reagents during execution of the method.
- Store the Neogen Molecular Detection Assay – *Listeria* Right Now as indicated on the package and in the product instructions.
- Always use the Neogen Molecular Detection Assay – *Listeria* Right Now by the expiration date.
- Always use the samplers contained in the kit. Do not use samplers from other sources.
- Right Now Sampler is single-use.
- Do not pause once reagent pellet rehydration has started. MDALRN96 does not allow the procedure to be interrupted at any step.
- It is strongly recommended that female laboratory staff be informed of the risk to a developing fetus resulting from infection of the mother through exposure to *Listeria monocytogenes*.

To reduce the risks associated with exposure to chemicals and biohazards:

- Perform pathogen testing in a properly equipped laboratory under the control of trained personnel.
- Users must clean food contact surface after sampling to remove any remaining swabbing solution.

To reduce the risks associated with cross-contamination while preparing the assay:

- Always wear gloves (to protect the user and prevent introduction of nucleases).

To reduce the risks associated with environmental contamination:

- Follow current industry standards for disposal of contaminated waste.
- Users must clean food contact surface after sampling to remove any residual swabbing solution.

**CAUTION**

To reduce the risks associated with exposure to hot liquids:

- Do not exceed the recommended temperature setting on the heater.
- Do not exceed the recommended heating time.
- Use an appropriate, calibrated thermometer to verify the Neogen® Molecular Detection Heat Block Insert temperature (e.g., a partial immersion thermometer or digital thermocouple thermometer, not a total immersion thermometer). The thermometer must be placed in the designated location in the Neogen Molecular Detection Heat Block Insert.

To reduce the risks associated with cross-contamination while preparing the assay:

- Change gloves prior to reagent pellet hydration.
- Use of long-tip, sterile, aerosol barrier (filtered), molecular biology grade pipette tips are recommended.
- Use a new pipette tip for each sample transfer.
- Use Good Laboratory Practices to transfer the sample to the lysis tube. To avoid cross-contamination between samples, the user may choose to clean the shaft of the pipette with a wipe with 10% bleach and then with a clean paper towel or use extended length pipette tips.
- Use a molecular biology workstation containing a germicidal lamp where available. Periodically decontaminate laboratory benches and equipment (pipettes, cap/decap tools, etc.) with a 1–5% (v:v in water) household bleach solution or RNA removal solution.

To reduce the risks associated with a false-positive result:

- Never open tubes post-amplification.
- Always dispose of the contaminated tubes by soaking in a 1–5% (v:v in water) household bleach solution for 1 hour, and away from the assay preparation area.
- Never autoclave reagent tubes post-amplification.

Consult the Safety Data Sheet for additional information and local, state, and national regulations for disposal.

If you have questions about specific applications or procedures, please visit our website at neogen.com or contact your local Neogen representative or distributor.

User Responsibility

Users are responsible for familiarizing themselves with product instructions and information. Visit our website at neogen.com or contact your local Neogen representative or distributor for more information.

When selecting a test method, it is important to recognize that external factors such as sampling methods, testing protocols, sample preparation, handling, laboratory technique, and the sample itself may influence results.

It is the user's responsibility in selecting any test method or product to evaluate a sufficient number of samples with the appropriate matrices and microbial challenges to satisfy the user that the chosen test method meets the user's criteria.

It is also the user's responsibility to determine that any test methods and results meet their customers' and suppliers' requirements.

As with any test method, results obtained from the use of any Neogen Food Safety product do not constitute a guarantee of the quality of the samples or processes tested.

To help customers evaluate the method for various food matrices, Neogen has developed the Neogen® Molecular Detection Assay – *Listeria* Right Now Matrix Control kit. When needed, use the Matrix Control (RNMC96) to determine that the matrix does not interfere with the Neogen Molecular Detection Assay – *Listeria* Right Now assay. Test several samples, representative of the environment matrix, i.e., samples obtained from different origins, during any validation period when adopting the Neogen method or when cleaning and sanitation protocols have changed.

Limitation of Warranties / Limited Remedy

EXCEPT AS EXPRESSLY STATED IN A LIMITED WARRANTY SECTION OF INDIVIDUAL PRODUCT PACKAGING, NEOGEN DISCLAIMS ALL EXPRESS AND IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO, ANY WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR USE. If any Neogen Food Safety Product is defective, Neogen or its authorized distributor will, at its option, replace or refund the purchase price of the product. These are your exclusive remedies. You must promptly notify Neogen within sixty days of discovery of any suspected defects in a product and return it to Neogen. Please contact your Neogen Food Safety representative or authorized Neogen distributor for any further questions.

Limitation of Neogen Liability

NEOGEN WILL NOT BE LIABLE FOR ANY LOSS OR DAMAGES, WHETHER DIRECT, INDIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES, INCLUDING BUT NOT LIMITED TO LOST PROFITS. In no event shall Neogen's liability under any legal theory exceed the purchase price of the product alleged to be defective.

Storage and Disposal

Store the Neogen® Molecular Detection Assay – *Listeria* Right Now at 2–8°C. Do not freeze. Keep the kit away from light during storage.

After opening the kit, check that the foil pouch, plastic, and glass bottles are undamaged. If the kit contents are damaged, do not use it.

After opening, unused reagent tubes should always be stored in the resealable pouch with the desiccant inside to maintain stability of the lyophilized reagents. Store resealed pouches at 2–8°C for no longer than 60 days.

After use, recap the Right Now™ Lysis Reagent Resuspension Buffer until firmly capped.

Reconstituted Right Now™ Lysis Buffer Reagent can be stored capped at 2–8°C for up to 21 days.

If necessary, the cap on the Lysis Buffer Reagent can be further secured with Parafilm®.

Do not use Neogen® Molecular Detection Assay – *Listeria* Right Now past the expiration date. The expiration date and lot number are noted on the outside label of the box.

When testing is complete, follow national, state, and local regulations for the disposal of contaminated waste. Consult the Safety Data Sheet for additional information.

Instructions for Use

Follow all instructions carefully. Failure to do so may lead to inaccurate results.

Periodically decontaminate laboratory benches and equipment (pipettes, cap/decap tools, etc.) with a 1–5% (v:v in water) household bleach solution or nucleic acid removal solution.

The user should complete the Neogen® Molecular Detection System operator qualification (OQ) training, as described in the “Installation Qualification (IQ) / Operational Qualification (OQ) Protocols and Instructions for Neogen Molecular Detection System” document (3).

1. Preparation of Right Now™ Sampler

- 1.1 Label the Right Now Sampler tube with the sample identification (ID).
- 1.2 Flick the Right Now Sampler tube twice to bring the sampling solution to the bottom of the tube.
- 1.3 Open the Right Now Sampler tube and squeeze the sampler to remove extra sampling solution.
- 1.4 Aseptically, pour the extra sampling solution out of the tube.
- 1.5 Return the Right Now Sampler to the empty tube.


2. Sample collection

- 2.1. Sample the surface by swabbing horizontally, vertically, and diagonally to provide full coverage of the test area (1×1 in to 12×12 in). Flip the Right Now Sampler over to the other side, and swab the same sampled surface.
- 2.2. Place the Right Now Sampler back into the empty tube.
- 2.3. Store the samplers at 2–8°C if not testing immediately. The Right Now Sampler can be held for up to 24 hours at 2–8°C before testing.


3. Preparation of Right Now™ Lysis Buffer Reagent (lytic enzyme)

- 3.1. Reconstitute one bottle of Right Now™ Lysis Buffer Reagent that contains the lyophilized lytic enzyme with 28 mL of the Right Now™ Lysis Reagent Resuspension Buffer. Swirl gently to mix. Allow the lysis buffer to warm up to room temperature before use.


NOTE: Rehydrated lysis buffer can be stored at 2–8°C for 21 days. Write the date on the vial that the lytic enzyme was prepared.

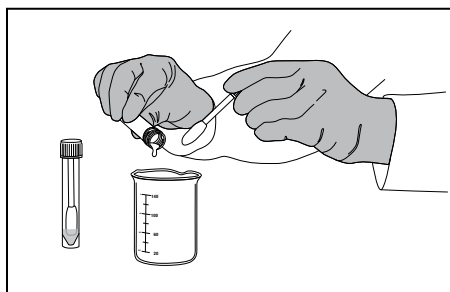Aseptically pour off the excess sampling solution in tube before swabbing.

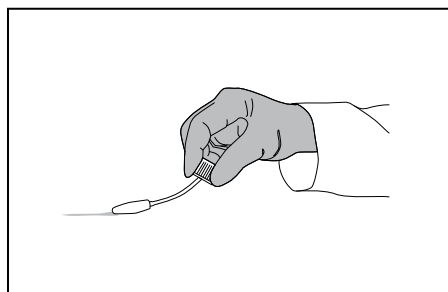Sample the surface by swabbing horizontally, vertically, and diagonally to provide full coverage of the test area (1×1 in to 12×12 in). Flip the Right Now Sampler over to the other side and swab the same sampled surface.

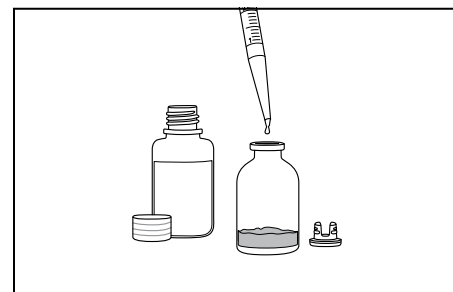Reconstitute one bottle of lysis buffer reagent (lyophilized lytic enzyme) with 28 mL of lysis reagent resuspension buffer. Swirl gently to mix. Allow the lysis buffer to warm up to room temperature before use.

Preparation of equipment

4. Preparation of the Neogen® Molecular Detection Speed Loader Tray

- 4.1 Wet a cloth or disposable towel with a 1–5% (v:v in water) household bleach solution and wipe the Neogen Molecular Detection Speed Loader Tray.
- 4.2 Rinse the Neogen Molecular Detection Speed Loader Tray with water.
- 4.3 Use a disposable towel to wipe the Neogen Molecular Detection Speed Loader Tray dry.
- 4.4 Ensure the Neogen Molecular Detection Speed Loader Tray is dry before use.

5. Preparation of the Neogen® Molecular Detection Chill Block Insert

- 5.1 Place the Neogen Molecular Detection Chill Block Insert directly on the laboratory bench. Use the block at ambient laboratory temperature.

6. Preparation of the Neogen® Molecular Detection Heat Block Insert

- 6.1 Place the Neogen Molecular Detection Heat Block Insert in a dry double block heater unit. Turn on the dry block heater unit and set the temperature to allow the Neogen Molecular Detection Heat Block Insert to reach and maintain a temperature of $100 \pm 1^\circ\text{C}$.

NOTE: Depending on the heater unit, allow approximately 30 minutes for the Neogen Molecular Detection Heat Block Insert to reach temperature. Using an appropriate, calibrated thermometer (e.g., a partial immersion thermometer, digital thermocouple thermometer, not a total immersion thermometer) placed in the designated location, verify that the Neogen Molecular Detection Heat Block Insert is at $100 \pm 1^\circ\text{C}$.

7. Preparation of the Neogen® Molecular Detection Instrument

- 7.1 Launch the Neogen® Molecular Detection Software and log in. Contact your Neogen Food Safety representative to ensure you have the most updated version of the software.
- 7.2 Turn on the Neogen Molecular Detection Instrument.
- 7.3 Create a New Run or edit a run by selecting a well for each sample. Refer to the Neogen Molecular Detection System User Manual for details.
- 7.4 For the *Listeria* Right Now (LRN) Assay, select the LRN icon on the Neogen Molecular Detection Software and assign to each well.
- 7.5 Assign the LRN NC icon and the LRN RC icon to corresponding wells within the run for the negative control and reagent control respectively.
- 7.6. LRN NC is set up for each run and one LRN RC is set up for the kit.

NOTE: The Neogen Molecular Detection Instrument must reach and maintain a temperature of 60°C before inserting the Neogen Molecular Detection Speed Loader Tray with reaction tubes. This heating step takes approximately 20 minutes and is indicated by an ORANGE light on the instrument's status bar. When the instrument is ready to start a run, the status bar will turn GREEN.


8. Lysis for *Listeria* Right Now™ Assay

- 8.1 Transfer 1 mL of the Right Now™ Lysis Buffer Reagent solution prepared in 3.1 into the Right Now Sampler used to collect the sample (section 2). Avoid touching the sides of the Right Now Sampler tube. Change tips with each sample to prevent contamination of the lysis buffer.


NOTE: Return the unused lysis buffer to the refrigerator after use (within 1 hour).

- 8.2 Close the tube cap. Mix the Right Now Sampler tube vigorously with a vortex mixer for 5–10 seconds to express the sample from the Right Now Sampler.
- 8.3 Place the number of individual cluster tubes in an empty rack. One empty cluster tube is required for each sample to be tested. Always include an additional empty cluster tube for a negative control (NC).
- 8.4 Remove the Right Now Sampler from the tube prepared in 8.1 and stand it on its cap on the bench. Avoid touching the sampler to prevent contamination.
- 8.5 Transfer 500 µL from each Right Now Sampler tube into a corresponding empty cluster tube using an extended length pipette tip. Return the sampler to the tube and secure the cap. Note: Store the Right Now Sampler tube with the remaining sample in the refrigerator for additional re-testing as needed.


NOTE: Extended-length pipette tips (that are longer than standard pipette tips) are recommended. In case the pipette shaft touches the Right Now Sampler tube, eject the pipette tip, then use a towel sprayed with 1–5% bleach to thoroughly wipe the shaft of the pipette to prevent cross-contact contamination between samples.

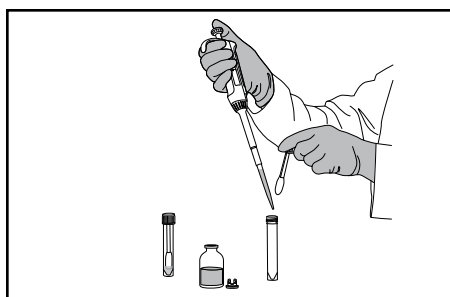Transfer 1 mL of the Right Now™ Lysis Buffer Reagent solution into the Right Now™ Sampler tube used to collect the sample. Avoid touching the sides of the Sampler tube.

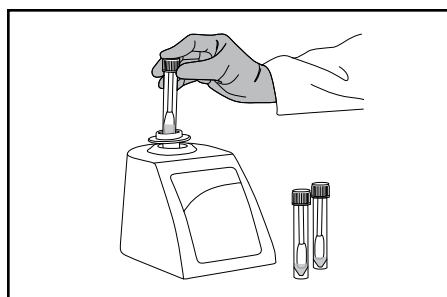
Close the tube cap. Mix the Sampler tube vigorously with a vortex mixer for 5–10 seconds to express the sample from the Sampler.

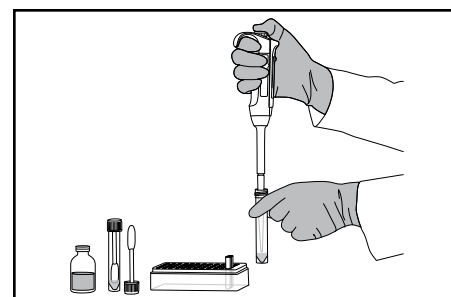
Using an extended-length pipette tip, transfer 500 µL of each sample from the Sampler tube into the corresponding empty cluster tube.

- 8.6 When all samples have been transferred, prepare the NC by transferring 500 µL of lysis buffer (prepared in 3.1) into an empty cluster tube. Do not use water as a NC.
- 8.7 Verify that the temperature of the Neogen Molecular Detection Heat Block Insert is at $100 \pm 1^\circ\text{C}$.
- 8.8 Place the uncovered rack of cluster tubes in the Neogen Molecular Detection Heat Block Insert and heat for 10 ± 1 minutes. During heating, the lysate may turn cloudy.
- 8.9 Remove the uncovered rack of cluster tubes from the heating block and allow it to cool in the Neogen Molecular Detection Chill Block Insert at least 5 minutes and a maximum of 10 minutes. The Neogen Molecular Chill Block Insert, used at ambient temperature, should sit directly on the laboratory bench. When cool, the lysate should be clear.

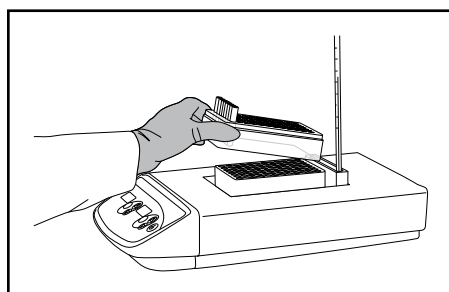
Place the uncovered rack of cluster tubes in the Neogen Molecular Detection Heat Block Insert and heat for 10 ± 1 minutes. During heating, the lysate may turn cloudy.

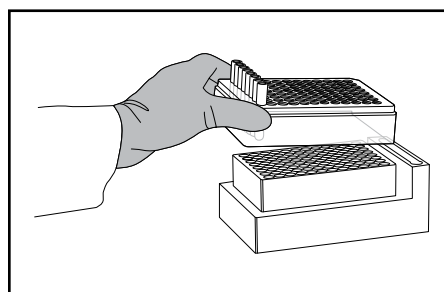
Remove the uncovered rack of cluster tubes from the heating block and allow it to cool in the Neogen Molecular Detection Chill Block Insert at least 5 minutes and a maximum of 10 minutes.

9. Amplification

- 9.1 One Neogen Molecular Detection Assay – *Listeria* Right Now Reagent Tube is required for each sample and the NC.
 - 9.1.1 Reagent tube strips can be cut to the desired tube number. Select the number of individual reagent tubes or 8-tube strips needed.
 - 9.1.2 Avoid disturbing the reagent pellets from the bottom of the tubes.
- 9.2 Remove one *Listeria* Right Now Reagent Control (RC) tube and place it in the rack.
- 9.3 For LRN NC lysate, select one *Listeria* Right Now Reagent Tube and place in the same rack. To avoid cross-contamination, decap one reagent tube strip at a time and use a new pipette tip for each transfer step.
- 9.4 Transfer the lysates to a *Listeria* Right Now Reagent Tube as described below.
- 9.5 Use the Neogen® Molecular Detection Cap/Decap Tool-reagent to decap the reagent tubes – one reagent tube strip at a time. Discard cap.
 - 9.5.1 Transfer 50 µL of sample lysate from the upper ½ of the liquid in the cluster tube (avoid precipitate) into the Neogen *Listeria* Right Now reagent tube. Dispense at an angle to avoid disturbing the pellets. Mix by gently pipetting up and down 5 times.**
 - 9.5.2 Repeat step 9.5.1 until all individual sample lysate has been added to a reagent tube in the strip.
 - 9.5.3 Cover the reagent tubes with the provided extra caps and use the rounded side of the Neogen Molecular Detection Cap/Decap Tool-reagent to apply pressure in a back-and-forth motion, ensuring that the cap is tightly applied.
 - 9.5.4 Repeat steps 9.5.1 to 9.5.3 as needed for the number of samples to be tested.
 - 9.5.5 When all sample lysates have been transferred, transfer 50 µL of the NC lysate to one *Listeria* Right Now Reagent Tube. Mix by gently pipetting up and down 5 times.
 - 9.5.6 Finally, transfer 50 µL of NC lysate into a Right Now™ Reagent Control (RC). Mix by gently pipetting up and down 5 times.

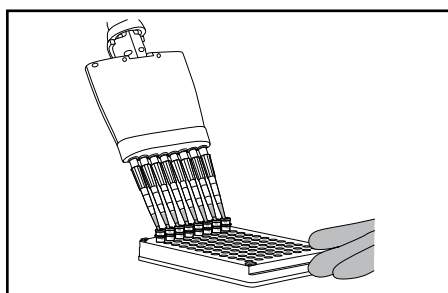Transfer 50 µL of sample lysate from the upper ½ of the cluster tube into the LRN Reagent Tube. Dispense at an angle to avoid disturbing the pellets. Mix by gently pipetting up and down 5 times.

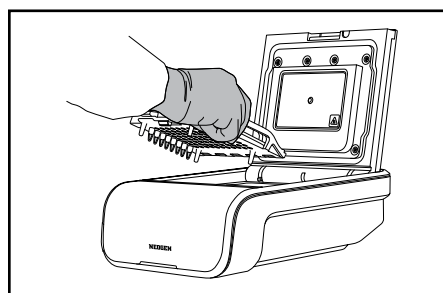Load capped tubes into a clean and decontaminated Neogen Molecular Detection Speed Loader Tray. Close and place it into the Neogen MDS. Close the lid to start the assay

- 9.6 Load capped tubes into a clean and decontaminated Neogen Molecular Detection Speed Loader Tray. Close and latch the Neogen Molecular Detection Speed Loader Tray lid.
- 9.7 Review and confirm the configured run in the Neogen Molecular Detection Software.
- 9.8 Click the start button in the software and select the instrument for use. The selected instrument's lid automatically opens.
- 9.9 Place the Neogen Molecular Detection Speed Loader Tray into the Neogen Molecular Detection Instrument and close the lid to start the assay. Results are provided within 45 minutes, although positives may be detected sooner.

9.10 After the assay is complete, remove the Neogen Molecular Detection Speed Loader Tray from the Neogen Molecular Detection Instrument and dispose of the tubes by soaking in a 1–5% (v:v in water) household bleach solution for 1 hour, and away from the assay preparation area.

NOTE: To minimize the risk of false positives due to cross-contamination, never open reagent tubes containing amplified DNA. Always dispose of sealed reagent tubes by soaking in a 1–5% (v:v in water) household bleach solution for 1 hour, and away from the assay preparation area.

Results and Interpretation

An algorithm interprets the light output curve resulting from the detection of the nucleic acid amplification. Results are analyzed automatically by the software and are color-coded based on the result. A Positive or Negative result is determined by analysis of a number of unique curve parameters. Positive results are reported in real-time, while Negative results will be displayed after the run is completed.

NOTE: Even a negative sample will not give a zero reading as the system and Neogen Molecular Detection Assay – *Listeria* Right Now amplification reagents have a “background” relative light unit (RLU) reading.

Table 2: Result symbols for samples and interpretation as provided by the Neogen Molecular Detection System Software.

Well Type	Well Result Symbol	Result	Interpretation
Sample		Positive	The sample is positive for <i>Listeria</i> spp.
Sample		Negative	The sample is negative for <i>Listeria</i> spp.
Sample		Inhibited	The sample matrix was inhibitory for the assay. A re-test may be required.
Sample		Error	No bioluminescence was detected. A re-test may be required.

References

1. NF EN ISO/IEC 17025. General requirements for the competence of testing and calibration laboratories.
2. NF EN ISO 7218. Microbiology of food and animal feeding stuffs – General requirements and guidance for microbiological examinations.
3. Neogen Installation Qualification (IQ) / Operational Qualification (OQ) Protocols and Instructions for Neogen Molecular Detection System.

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© Neogen Corporation 2025. All rights reserved. Neogen is a registered trademark of Neogen Corporation.
All other brand and product names are trademarks or registered trademarks of their respective companies. Patent: neogen.com/patents
Neogen Corporation
Neogen Europe Ltd.
Neogen Ireland, Ltd.
Neogen Corporation, 620 Lesher Place, Lansing, MI 48912 USA.