-  (EN) Hazelnut Protein Rapid Kit
-  (FR) Kit Rapide Protéine Noisette
-  (DE) Haselnuss Protein Rapid Kit
-  (IT) Kit rapido per la rilevazione delle proteine di nocciola
-  (ES) Kit rápido para proteína de avellana
-  (NL) Hazelnoot Proteïne Sneltest
-  (SV) Hazelnut Protein Rapid Kit
-  (DA) Hasselnød Protein Rapid Kit
-  (NO) Hasselnøttprotein hurtigkit
-  (PT) Kit Rápido para Proteína de Avelã
-  (EL) Κιτ Ταχείας Δοκιμασίας Πρωτεΐνης Φουντουκιού
-  (PL) Szybki test na wykrywanie białek obecnych w orzechach laskowych
-  (RU) Rapid набор (протеин фундука)
-  (TR) Fındık Proteini Hızlı Test Kiti
-  (JA) ヘーゼルナッツプロテイン迅速測定キット
-  (ZH) 榛子蛋白快速检测试剂盒
-  (TH) ชุดทดสอบโปรตีนของเฮเซลนัทแบบรวดเร็ว
-  (KO) 헤이즐넛 단백질 신속 키트
-  (ID) Rapid Kit Protein Hazelnut
-  (AR) طقم سريع بروتين البندق

Product Instructions

Hazelnut Protein Rapid Kit

Lateral Flow Devices (LFD) for qualitative analysis of hazelnut proteins.

Product Description and Intended Use

The Neogen® Hazelnut Protein Rapid Kit is intended for screening for the presence of hazelnut proteins in clean-in-place (CIP) final rinse water, environmental swab samples, food ingredients and processed food products.

The Neogen Hazelnut Protein Rapid Kit utilizes a lateral flow device (LFD) that is an immunochromatographic test method utilizing antibodies which are specific for the detection of hazelnut proteins. Positive results are visualized by the presence of three lines: a control line, a hook line and a test line, when the hazelnut protein is present in the range of concentrations from a lower level of 2 ppm and a higher level of approximately 10 mg/mL of hazelnut protein in saline solution (0.9% sodium chloride). These limits may vary depending on the matrix tested. Samples containing more than 10% hazelnut may result in an invalid test which will result in the absence of the hook line and/or test line.

The Neogen Hazelnut Protein Rapid Kit is intended for use in the food and beverage industry by trained personnel. Neogen has not documented the use of this product in industries other than food or beverage. For example, Neogen has not documented this product for testing pharmaceutical, cosmetic, clinical or veterinary samples. The Neogen Hazelnut Protein Rapid Kit has not been evaluated with all possible food products, food processes and testing protocols.

Neogen Hazelnut Protein Rapid Kit contains 25 tests, described in Table 1.

Table 1. Kit components

Item	Identification	Quantity	Storage
Neogen® Hazelnut Protein Lateral Flow Device (LFD)	Lateral flow device in a plastic cassette	25 devices individually packed.	Store at 2-8°C. Do not freeze.
Neogen® Extraction Buffer	Bottle with Extraction Buffer	1 bottle with 50 mL	Store at 2-8°C. Do not freeze.
Dilution Tubes	Microcentrifuge tube (2.2 mL volume capacity)	26 tubes	Store in a clean dry place.

Materials not provided in the kit:

- Swabs and pipettes.
- The use of vortex, timer and balance are recommended but not required for all samples.
- The use of a centrifuge is required for all chocolate and gum samples and it is recommended, but not required for all solid samples.

Safety

The user should read, understand, and follow all safety information in the instructions for the Neogen Hazelnut Protein Rapid Kit. Retain the safety instructions for future reference.

⚠ WARNING: Indicates a hazardous situation, which, if not avoided, could result in death or serious injury and/or property damage.

NOTICE: Indicates a potentially hazardous situation, which, if not avoided, could result in property damage.

⚠ WARNING

To reduce the risks associated with inaccurate results:

- Neogen has not documented the use of Neogen Hazelnut Protein Rapid Kit in industries other than food or beverage. For example, Neogen has not documented this product for testing pharmaceutical, cosmetic, clinical or veterinary samples.
- Neogen Hazelnut Protein LFD should be read 11 ± 1 minutes after sample has been loaded on the lateral flow device.
- The Neogen Extraction Buffer is designated for use with a specific lot of Neogen Hazelnut Protein LFD. Do NOT interchange Neogen Hazelnut Protein Rapid Kit components with other lots or kits.
- The Neogen Extraction Buffer is designated for use with a specific lot of Neogen Hazelnut Protein LFD. Dispose of any remaining Neogen Extraction Buffer once all Neogen Hazelnut Protein Lateral Flow Devices have been used.

- Store the Neogen Hazelnut Protein Rapid Kit as indicated on the package and in the product instructions.
- Always use the Neogen Hazelnut Protein Rapid Kit by the expiration date.
- Always use the Neogen Hazelnut Protein Rapid Kit at 20-25°C temperature.
- Neogen® Allergen Protein Testing Kits are not intended for the detection of hydrolyzed proteins.
- The Neogen Allergen Protein Testing Kits are designed to detect proteins from processed food once they are solubilized into Neogen Extraction Buffer. Some food processing methods may limit detection of these target proteins.
- Some food processing may affect the detection of food proteins with Neogen Allergen Protein Testing Kits. Users should verify that the method is fit for purpose to meet user's requirements.

To reduce the risks associated to false negative results:

- Use the Neogen Hazelnut Protein Rapid Kit for food and environmental samples that have been validated internally or by a third party.

To reduce the risks associated with exposure to chemicals:

- Neogen Hazelnut Protein Rapid Kit is intended for use in the food and beverage industries by trained personnel.

NOTICE**To reduce the risks of inaccurate results:**

- Refer to the Interpretation of Results section of the product instructions, to ensure accurate interpretation of the Neogen Hazelnut Protein LFD.

Consult the Safety Data Sheet for additional information.

For information on documentation of product performance, visit our website at www.neogen.com or contact your local Neogen representative or distributor.

User Responsibility

Users are responsible for familiarizing themselves with product instructions and information. Visit our website at www.neogen.com, or contact your local Neogen representative or distributor for more information.

As with all test methods used for food analysis the test matrix can influence the results. When selecting a test method, it is important to recognize that external factors such as sampling methods, testing protocols, sample preparation, handling, and laboratory technique may influence results. The food sample itself may influence results.

It is the user's responsibility in selecting any test method or product to evaluate a sufficient number of samples to satisfy the user that the chosen test method meets the user's criteria.

It is also the user's responsibility to determine that any test methods and results meet its customers' and suppliers' requirements.

As with any test method, results obtained from use of any Neogen Food Safety product do not constitute a guarantee of the quality of the matrices or processes tested.

Limitation of Warranties/Limited Remedy

EXCEPT AS EXPRESSLY STATED IN A LIMITED WARRANTY SECTION OF INDIVIDUAL PRODUCT PACKAGING, NEOGEN DISCLAIMS ALL EXPRESS AND IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO, ANY WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR USE. If any Neogen Food Safety Product is defective, Neogen or its authorized distributor will, at its option, replace or refund the purchase price of the product. These are your exclusive remedies. Please contact your Neogen representative or authorized Neogen distributor for any further questions.

Limitation of Neogen Liability

NEOGEN WILL NOT BE LIABLE FOR ANY LOSS OR DAMAGES, WHETHER DIRECT, INDIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES, INCLUDING BUT NOT LIMITED TO LOST PROFITS. In no event shall Neogen's liability under any legal theory exceed the purchase price of the product alleged to be defective.

Storage and Disposal

Store all Neogen Hazelnut Protein Rapid Kit components at 2-8°C.

Neogen Hazelnut Protein Rapid Kit components should not be frozen, exposed to UV light or exposed to prolonged heat (>30°C).

Neogen Hazelnut Protein Rapid Kit components should not be used past the expiration date. Expiration date and lot number are noted on the outside label of the box.

Please note that each Neogen Extraction Buffer lot is validated specifically for each LFD lot and is not interchangeable with any other lots or kits.

Dispose according to current local/regional/industry standards and regulations.

Instructions for Use

Follow all instructions carefully. Failure to do so may lead to inaccurate results. Ensure that all Neogen Hazelnut Protein Rapid Kit components are at ambient temperature (20-25°C) prior to use.


Sample Analysis

1. CIP Final Rinse Water Samples

- 1.1 Label one microcentrifuge tube for each CIP sample.
- 1.2 Add 800 μ L of Neogen Extraction Buffer to a labeled microcentrifuge tube.
- 1.3 Add 200 μ L of CIP final rinse water sample. Shake vigorously or vortex for 15 seconds to mix thoroughly to obtain an extracted sample.


Note: The pH of extracted sample should be between 5 and 10. Proceed to Troubleshooting section for additional information.

- 1.4 Remove one Neogen Hazelnut Protein LFD from package and place on a clean, dry, flat surface.
- 1.5 Transfer 100 μ L of the extracted sample prepared in 1.3 using a clean pipette or pipette tip and apply it to the sample well on the Neogen Hazelnut Protein LFD. Start the timer for 11 \pm 1 minutes. Proceed to Interpretation of Results section.

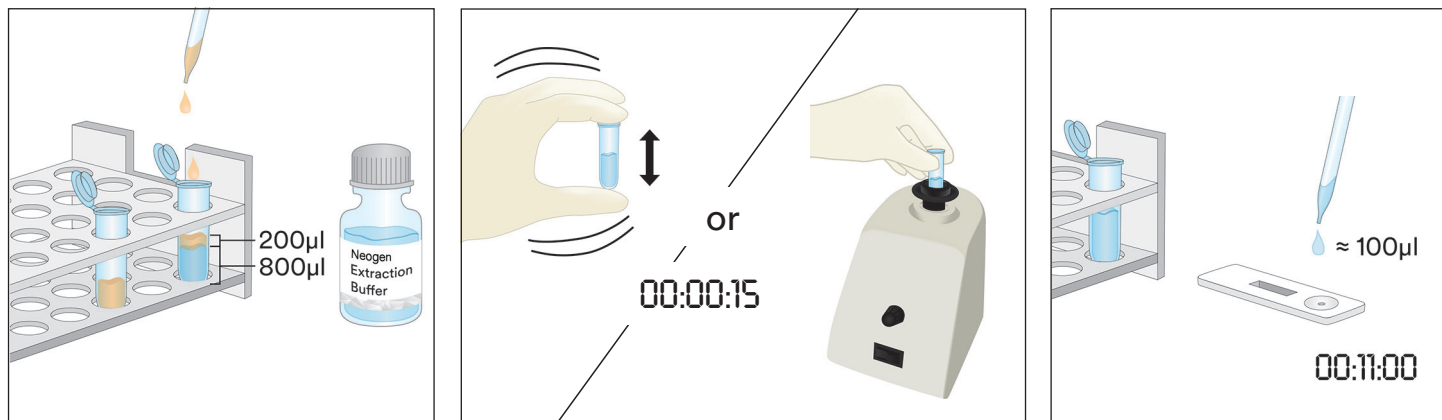2. Environmental Swab Samples

- 2.1 Label one microcentrifuge tube for each Environmental Swab Sample.
- 2.2 Add 500 μ L of Neogen Extraction Buffer into a labeled microcentrifuge tube.
- 2.3 Take a clean swab and dip the entire tip into the microcentrifuge tube wetting the tip with Neogen Extraction Buffer. Gently express excess liquid from tip by pressing swab tip lightly on the inside of the tube.
- 2.4 Take wetted swab and survey a 10 X 10 cm surface area maintaining the swab at a 30° angle with the surface. Rub the swab slowly and thoroughly over the surface area. Rub the swab three times over this surface, reversing direction between alternating strokes.
- 2.5 Take swab and insert it back into the pre-labeled tube and swirl swab several times to release any residues that might be on the surface of swab into the Neogen Extraction Buffer. Break off swab tip in tube, cap tightly and mix well to obtain an extracted sample.
- 2.6 Remove one Neogen Hazelnut Protein LFD from the package and place on a clean, dry, flat surface.
- 2.7 Transfer 100 μ L of the extracted sample prepared in 2.5 using a clean pipette or pipette tip and apply it to the sample well on the Neogen Hazelnut Protein LFD. Start the timer for 11 \pm 1 minutes. Proceed to Interpretation of Results section.

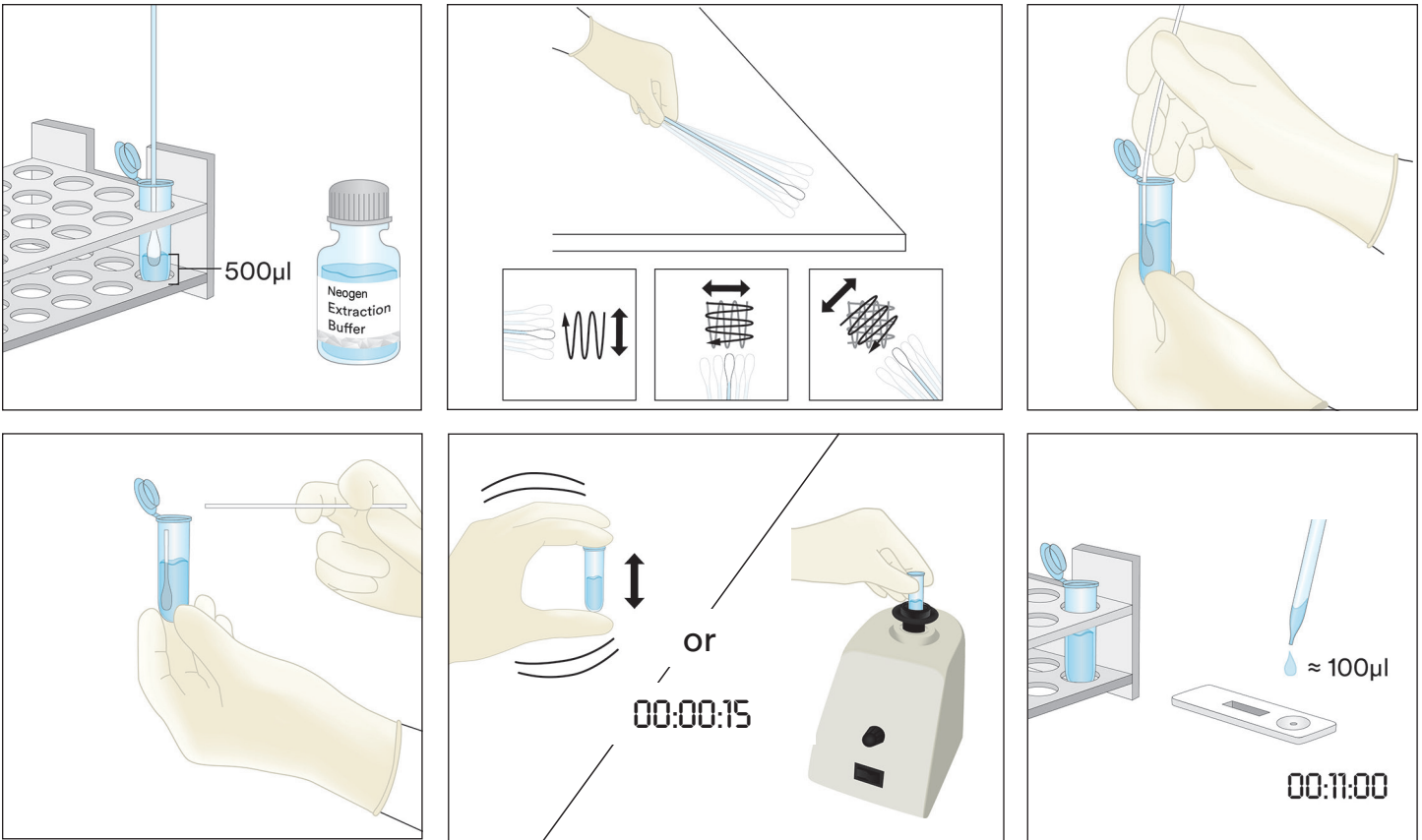
3. Liquid Samples Except Liquid Chocolate Samples

- 3.1 Label one microcentrifuge tube for each Liquid sample.
- 3.2 Measure 900 µL of Neogen Extraction Buffer into a labeled microcentrifuge tube.
- 3.3 Add 100 µL of a well-mixed sample. Shake vigorously or vortex for 15 seconds to mix thoroughly to obtain an extracted sample.


Note: The pH of extracted sample should be between 5 and 10. Proceed to Troubleshooting section for additional information.

- 3.4 Remove one Neogen Hazelnut Protein LFD from package and place on a clean, dry, flat surface.
- 3.5 Transfer 100 µL of the extracted sample prepared in 3.3 from the middle (aqueous) layer using a clean pipette or pipette tip and apply it to the sample well on the Neogen Hazelnut Protein LFD. Start the timer for 11 ± 1 minutes. Proceed to Interpretation of Results section.

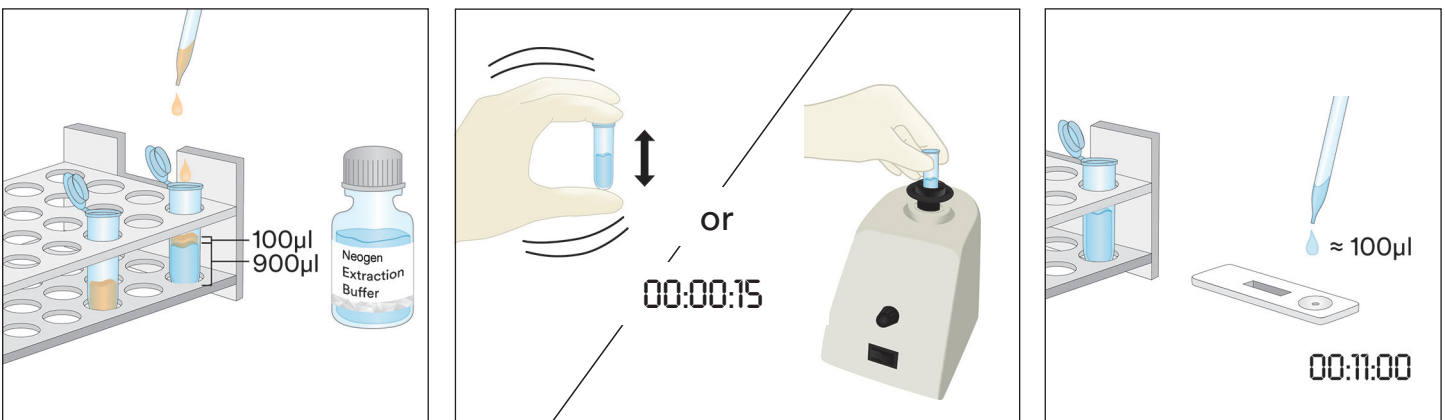4. Liquid Chocolate Samples

- 4.1 Label two microcentrifuge tubes for each Liquid Chocolate Sample.
- 4.2 To make a prepared sample, add 500 μL of a well-mixed Liquid Chocolate Sample and add 500 μL of Neogen Extraction Buffer into one labeled microcentrifuge tube and shake vigorously to mix thoroughly or vortex for approximately 15 seconds.
- 4.3 Add 900 μL of Neogen Extraction Buffer into the second labeled microcentrifuge tube and add 100 μL of the Prepared Sample from step 4.2. Shake vigorously to mix thoroughly or vortex for approximately 15 seconds.
- 4.4 Centrifuge for 20-30 seconds at 5000-7000 rpm (3000 x g). The supernatant is the extracted sample.
- 4.5 Remove one Neogen Hazelnut Protein LFD from package and place on a clean, dry, flat surface.
- 4.6 Transfer 100 μL of the extracted sample prepared in 4.4 from the middle (aqueous) layer using a clean pipette or pipette tip and apply it to the sample well on the Neogen Hazelnut Protein LFD. Start the timer for 11 ± 1 minutes. Proceed to Interpretation of Results section.

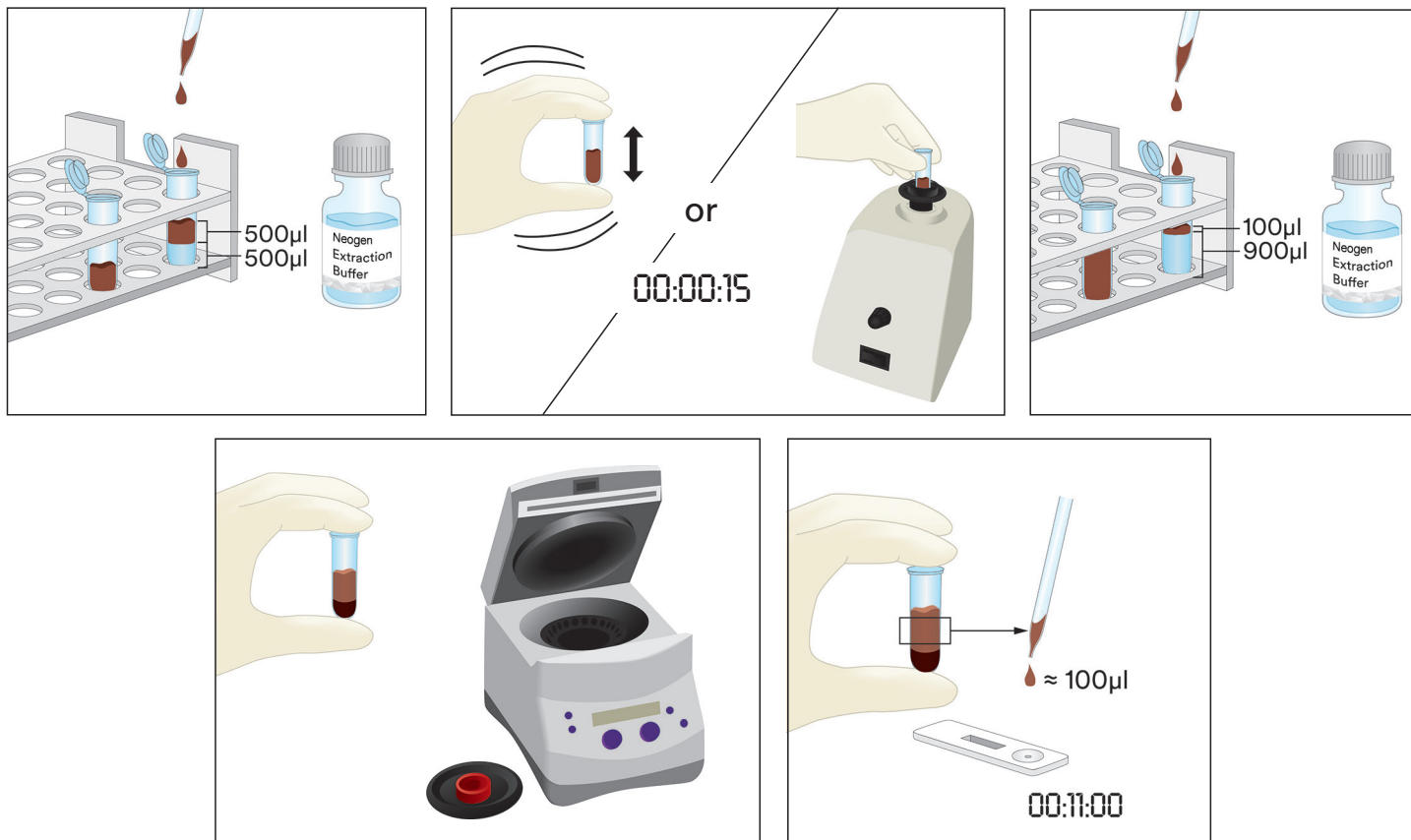5. Solid Samples

- 5.1 Label one microcentrifuge tube for each non-liquid sample.
- 5.2 Grind sample into a fine homogeneous powder.
- 5.3 Measure 0.2 g of sample into a labeled microcentrifuge tube.
- 5.4 Add 1.8 mL of the Neogen Extraction Buffer to the sample in the microcentrifuge tube. Shake vigorously or vortex for 15 seconds to mix thoroughly to obtain an extracted sample.
Note: Melt chocolate samples.
- 5.5 Centrifuge for 20-30 seconds at 5000-7000 rpm (3000 x g). The supernatant is the extracted sample.
- 5.6 Remove one Neogen Hazelnut Protein LFD from package and place on a clean, dry, flat surface.
- 5.7 Transfer 100 μL of the extracted sample prepared in 5.5 from the middle (aqueous) layer using a clean pipette or pipette tip and apply it to the sample well on the Neogen Hazelnut Protein LFD. Start the timer for 11 ± 1 minutes. Proceed to Interpretation of Results section.


Interpretation of Results

The control line is next to the letter C on the Neogen Hazelnut Protein LFD. The test line is next to the letter T on the Neogen Hazelnut Protein LFD. The hook line is next to the letter H and/or in between the T and C lines. Note: Some LFD cassettes might not have an H marking on the side of the plastic cassette.

Read LFD at 11 ± 1 minutes after application of the sample. A sample is considered to be:

- Negative for hazelnut protein when only the two lines furthest from the well, the control and the hook line, are visible on the Neogen Hazelnut Protein LFD.
- Positive for hazelnut protein when all three lines (test, hook and control line) are visible on Neogen Hazelnut Protein LFD.
- Invalid, if the Neogen Hazelnut Protein LFD does not develop the middle hook line. The sample may be overloaded with hazelnut protein and needs to be diluted.
- Invalid, if the Neogen Hazelnut Protein LFD does not develop the control line.

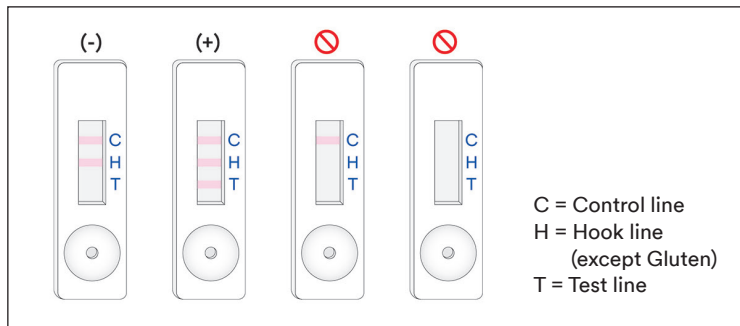Any reading after 12 minutes from the initial application of the sample into the Neogen Hazelnut Protein LFD should be considered invalid. A reading at this time cannot be interpreted and can lead to erroneous results.

Troubleshooting

1. Sample fails to migrate across the strip within the first 5 minutes after application of the sample in the Neogen Hazelnut Protein LFD.
The sample may be too viscous and needs to be centrifuged if this was not already done during the preparation of the sample. If sample was already centrifuged, then preparing a 1:1 dilution with the Neogen Extraction Buffer may be necessary. (Note: This may reduce the sensitivity to ~6 ppm for some matrices.)
2. A red dot appears on the test line but remainder of test line does not change color.
Sample particulate may have passed around the filter in the cassette, simply re-run the sample by taking a new Neogen Hazelnut Protein LFD from the kit and repeat the test.
3. The pH of extracted sample should be between 5 and 10. If pH is outside this range, further dilution may be required (i.e., prepare a 1:1 dilution with 100 μ L of the extracted sample and 100 μ L of Neogen Extraction Buffer. This may reduce the sensitivity to ~6 ppm for some matrices.)

If you have questions about specific applications or procedures, please contact your Neogen Food Safety representative or distributor.

Minimum Performance Characteristics

Lowest Limit of Detection ^(a)	2 ppm
Upper Limit of Detection	~10 mg/mL

- ^(a) The lowest limit of detection is defined as the lowest concentration of the allergen in a test sample that can be distinguished from a true blank sample at a specified probability level¹.

References

1. Abbott, M., Hayward, S., Ross, W., Godefroy, S.B., Ulberth, F., Van Hengel, A. J., Roberts, J., Akiyama, H., Popping, B., Yeung, J.M., Wehling, P., Taylor, S., Poms, R.E., and Delahaut, P. (2010). Appendix M: Validation Procedures for Quantitative Food Allergen ELISA Methods: Community Guidance and Best Practices. *J. AOAC Int.* 93, 442-450.

Explanation of Symbols

info.neogen.com/symbols

Neogen Food Safety

Neogen Corporation

620 Leshar Place
Lansing, MI 48912 USA
Neogen.com

Neogen Europe Ltd.

The Dairy School
Auchincruive
Ayr, KA6 5HU
Scotland, UK

Neogen Ireland, Ltd.

Bray Business Park, Bray
Co. Wicklow
A98YV29, Ireland

Neogen Corporation

620 Leshar Place Lansing, MI 48912 USA
www.neogen.com

Neogen is a trademark of Neogen Corporation.
© Neogen Corporation 2024.
Unauthorized use prohibited. All rights reserved.

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