

A young girl with long, wavy blonde hair is looking down at a bowl of fresh raspberries. She has a slight smile and is wearing a white top. The background is softly blurred, showing a warm, indoor setting.

Surface hygiene testing made simple

Neogen Clean-Trace™ Surface Protein Swabs

Monitoring the effectiveness of cleaning in food production and preparation areas can help you to reduce the risk of a food safety incident. Clean-Trace Surface Protein Swabs are a simple, quick and reliable method of checking for the removal of protein residues after cleaning.

Reliable results consistently and quickly

Clean-Trace™ Surface Protein Plus

Designed to detect residues of protein and reducing agents on surfaces using the biuret reaction.

- Semi-quantitative protein test with swab-click-read simplicity
- Easily access hard-to-reach areas
- Results available in 10 minutes or less if the surface contamination is high
- No instrument required
- Based on colour change for ease of interpretation
- Detects as little as 50µg of protein

Clean-Trace™ Surface Protein (Allergen) Test

Designed to detect residues of protein and reducing agents on surfaces and in process water using the biuret reaction to indicate the presence of potential allergens.

- Highly sensitive: 3 µg of proteins
- Fast: testing takes only 15 minutes
- Swab, click and incubate (55°C, 15 minutes), read
- Storage at ambient temperature
- Colour change makes it easy to read results

Tested Allergen	Minimal Detection
Egg	10-20 µg
Peanut	25-50 µg
Milk	10-20 µg
Gluten	40 µg
Soy	20-40 µg
Almond	20-40 µg




Easy to use

1 Swab

Swab surface firmly. If the surface is dry, wet the swab with the moisturiser supplied in the kit.


2 Click

Activate the swab and shake from side to side for 5 seconds.

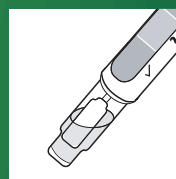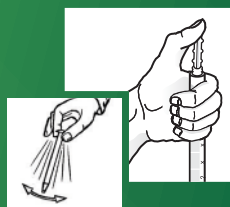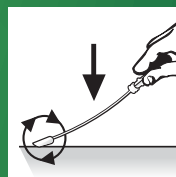

Protein Plus 3 Read

Leave the swab upright at room temperature for 10 minutes. Read the result.

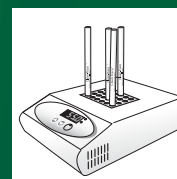

Allergen Test 3 Heat

Place the swab into a Digital Heating Block (set to 55°C) for 15 minutes.

4 Interpret results

Clean

Dirty

Product Code	Description	No. of tests
PRO100	Clean-Trace Surface Protein Plus	100
MPRO50	Clean-Trace Surface Protein Plus	50
ALLTEC60	Clean-Trace Surface Protein (Allergens)	60



To learn more about how Neogen products can help to protect your brand and business contact your local Neogen representative or visit neogenaustralasia.com.au